

Application News

No. C158

Direct Probe Ionization Mass Spectrometer

15 Second Screening Analysis of Cyanide in Blood Serum Without Pretreatment

Recent years have witnessed an increasing trend in incidents of crime and poisoning involving various legal drugs and toxic substances. The diversity in the types of used substances has led to such incidents becoming a social problem. In the fields of forensic medicine, forensic toxicology, and critical care, finding and identifying causative agents is a problem that requires establishment of an analysis method that is both quick and highly sensitive. There is also growing demand in these workplaces to further simplify the complex pretreatment processes and instrument operations as well as to increase analysis speed. While various analysis instruments have been utilized until now for analyzing specific components in blood, most instruments require complex pretreatment, such as extracting the target component from blood. What is needed is a screening method that best reduces the time and labor required to perform analysis.

Probe electrospray ionization (PESI) is a direct ionization technique that ionizes sampled target components by sampling samples using an ultrafine and minimally invasive probe and applying high voltage to the probe tip. This technique enables sample analysis without the need for a chromatograph (Fig. 1).

The DPiMS-2020, which combines PESI with a mass spectrometer, is suitable for simple screening analysis because it enables quick analysis of target components in samples without pretreatment regardless of whether samples are in liquid or solid form.

This article introduces a rapid screening method for detecting cyanide in blood serum that does not require pretreatment by utilizing the DPiMS-2020 and In-Source CID.

T. Murata

■ Analysis Sequence Without Pretreatment

Potassium cyanide was added to blood serum to obtain a final concentration of 10 µg/mL and then taurine and naphthalene dialdehyde were added to perform derivatization.^{*1} The obtained cyanide derivative (Fig. 2) was added to blood serum and used as the sample.

While complex pretreatment, such as that shown in Fig. 3, is required in conventional blood serum analysis, analysis that utilizes PESI can be performed using blood serum that contains cyanide derivatives either as-is or diluted with water by injecting it onto a small (10 µL) sample plate and setting the sample plate in the instrument.

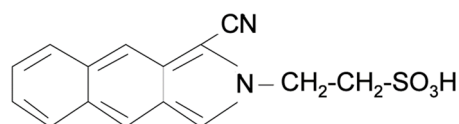


Fig. 2 Cyanide Derivative (MW 300)

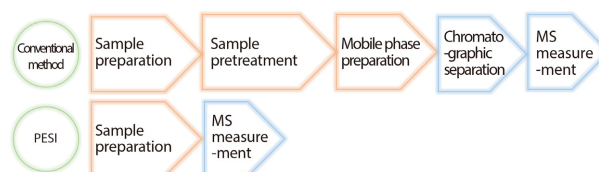


Fig. 3 Analysis Sequence Without Pretreatment

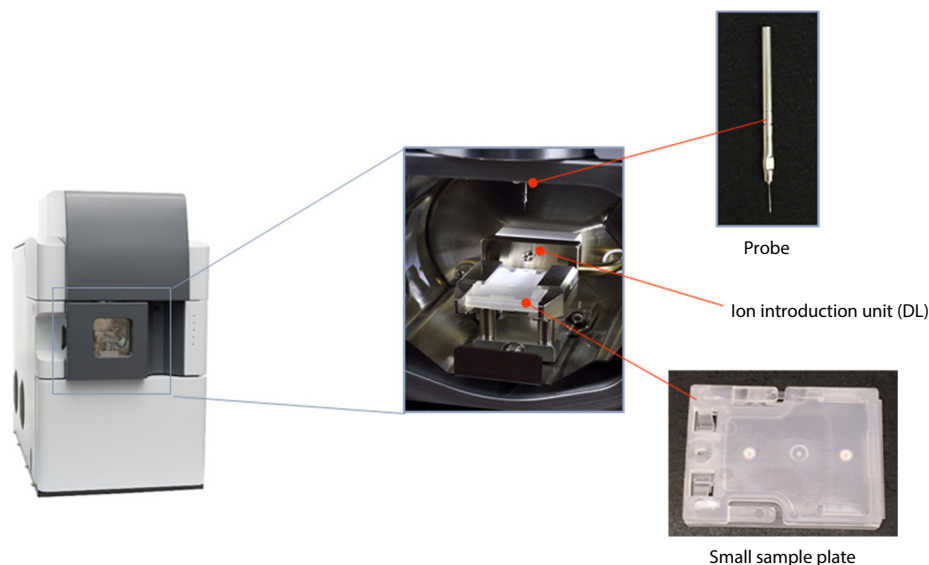


Fig. 1 DPiMS-2020

Structural Analysis Using In-Source CID Analysis

While a triple quadrupole mass spectrometer is usually used for structural analysis, a single quadrupole mass spectrometer can also obtain molecular structure information as well as molecular weight information by setting the lens system to a high voltage.

Analysis that utilizes PESI results in a unique mass chromatogram (Fig. 4) because the probe is driven at a constant frequency to repeat a process of sampling followed by ionization by applying a high voltage.

Fig. 5 shows the mass chromatogram obtained in our example. Applying a voltage of -80 V to DL bias and Q-array bias allows molecular structure information to be obtained and enables quick and simple screening analysis for cyanide in blood serum.

For reference, Fig. 6 shows the product ion (MS/MS) mass spectrum of the cyanide derivative obtained using the LCMS-8040 triple quadrupole mass spectrometer.

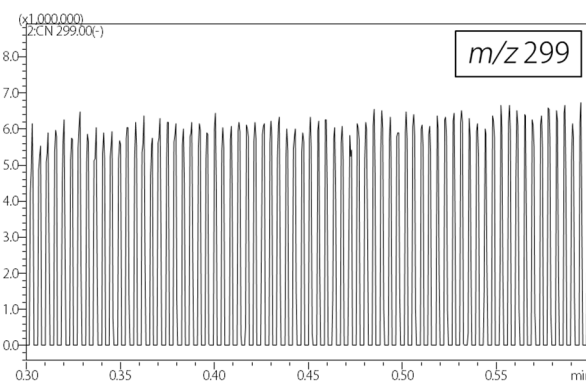


Fig. 4 Mass Chromatogram from DPiMS-2020

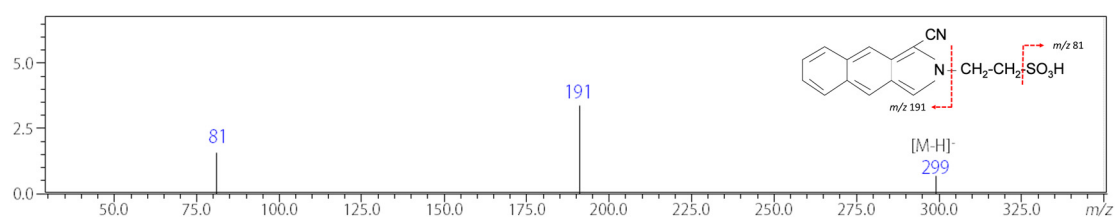


Fig. 5 Mass Spectrum from DPiMS-2020

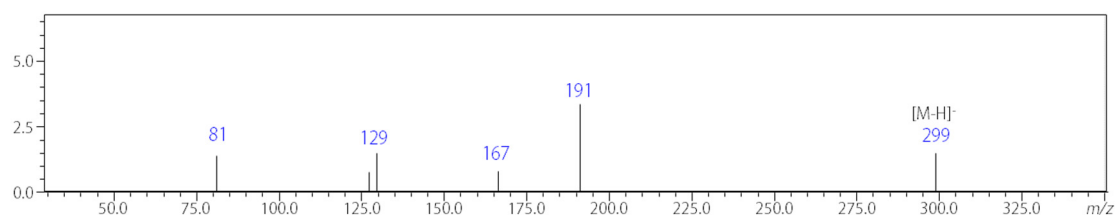


Fig. 6 Product Ion Mass Spectrum from LCMS-8040 (Reference)

DPiMS-2020 Analysis Conditions

When performing analysis using the DPiMS-2020, the drive conditions of the PESI probe and the analysis conditions of the mass spectrometer must be set.

Table 1 and 2 list the drive and analysis conditions respectively.

Table 1 PESI Drive Conditions

Ionization position	: -37 mm
Ionization stop time	: 100 msec
Sampling position	: -46 mm
Sampling stop time	: 50 msec
Probe speed	: 250 mm/s
Probe acceleration	: 0.63 G

Table 2 Mass Spectrometer Analysis Conditions

DL temperature	: $250\text{ }^{\circ}\text{C}$
Heater block temperature	: $35\text{ }^{\circ}\text{C}$
Interface voltage	: -2.45 kV (ESI – Negative mode)
DL bias voltage	: -80 V ($m/z\ 299$)
Q-array bias voltage	: -80 V ($m/z\ 299$)

References

*1 S. Chinaka, N. Takayama, Y. Michigami, and K. Ueda. *J. Chromatogr. B.* 713: 353–359 (1998)

Acknowledgments

We would like to thank associate professor Kei Zaitzu and assistant professor Yumi Hayashi at the Nagoya University Graduate School of Medicine for their guidance regarding data acquisition and sample preparation.

The product described in this document has not been approved or certified as a medical device under the Pharmaceutical and Medical Device Act of Japan. It cannot be used for the purpose of medical examination and treatment or related procedures.

First Edition: Jul. 2017



For Research Use Only. Not for use in diagnostic procedure.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Company names, product/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation or its affiliates, whether or not they are used with trademark symbol "TM" or "®". Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services. Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

Shimadzu Corporation

www.shimadzu.com/an/